

Fitting Instructions

Page 1

Make	Ford	Model	Tourneo Connect	Derivative	MK III 5 Seat	REV B
-------------	------	--------------	-----------------	-------------------	---------------	-------



Fitting Instructions						
Page 2						
Make	Ford	Model	Tourneo Connect	Derivative	MK III 5 Seat	REV B



Pre-installation.



Remove the tie down loops from the boot sill (T30 torx).



Unclip and remove the sill plastics.



Remove the two screws shown in the side panel (t30 torx).



Unclip and remove the side panel. The carpet can be folded out of the way.

Fitting Instructions					
Page 3					
Make	Ford	Model	Tourneo Connect	Derivative	MK III 5 Seat REV B



The upper, outer and inner legs are attached to the a-frame as shown.



Outer Leg:
4x M10x20mm Countersunk Set Screw



Inner Leg:
2x M10x20mm Countersunk Set Screw



Upper Leg:
4x M10x20mm Countersunk Set Screw



Temporarily attach the a-frame to the 2x threaded holes shown using the 2x M6x30mm Countersunk Set Screws.

Fitting Instructions						
Page 4						
Make	Ford	Model	Tourneo Connect	Derivative	MK III 5 Seat	REV B



Use the upper leg fixing hole as a marker to drill a 10mm clearance hole through the d-pillar.
Ensure care is taken to not drill too far through.



Attach the threaded pylon using the backplate shown.



Ensure the use of the washer.



The upper leg will attach to the vehicle as shown but do not attach at this point as the a-frame needs removing later on.




Before drilling for the inner and outer legs, check whether there is a wiring loom under the vehicle. This has only been seen on PHEV models so far. If there, it can be unclipped to allow for holes to be drilled and then re-attached.

Fitting Instructions						
Page 5						
Make	Ford	Model	Tourneo Connect	Derivative	MK III 5 Seat	REV B



Whilst the a-frame is in position, drill 8mm clearance holes through the floor for the inner and outer leg attachments. Then the a-frame can be removed in order to cut and re-fit the plastics.

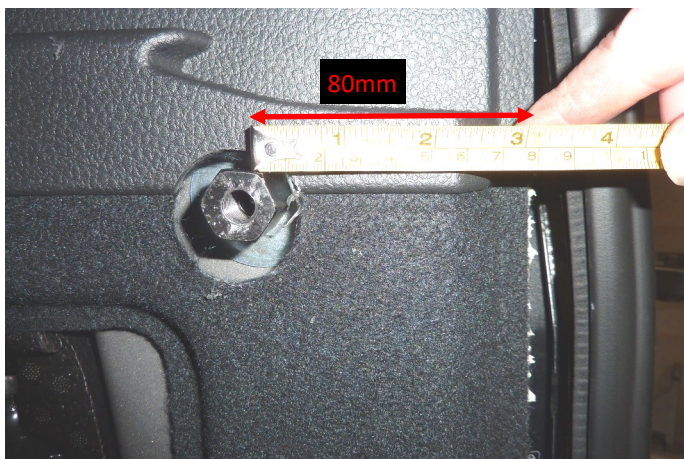


Cut the necessary hole in the side plastics for the pylon to protrude.



The photos and measurements can be used as an approximate guide.



Replace the side panel and use a pylon flange to make good the cut hole.

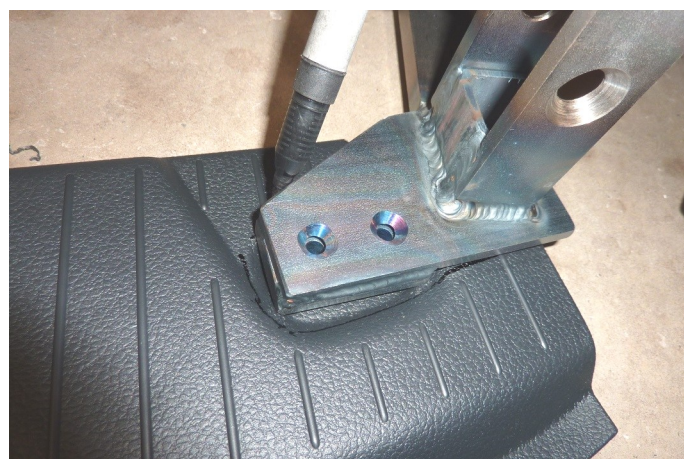

The outer leg can be used to draw around as a template for the necessary cut out in the plastic boot sill. The carpet should not need to be cut, it can fit over the leg brackets.

Fitting Instructions						
Page 6						
Make	Ford	Model	Tourneo Connect	Derivative	MK III 5 Seat	REV B



Fit the a-frame to the vehicle. A small cut out will be required for the inner leg as shown. This can be cut using a utility knife in situ.



The inner leg uses a spreader plate under the vehicle. Apply an appropriate corrosion resistant sealant.



Fitting Instructions

Page 7

Make	Ford	Model	Tourneo Connect	Derivative	MK III 5 Seat	REV B
------	------	-------	-----------------	------------	---------------	-------



Run the power cable to the 12V battery which is at the left hand side of the engine bay.

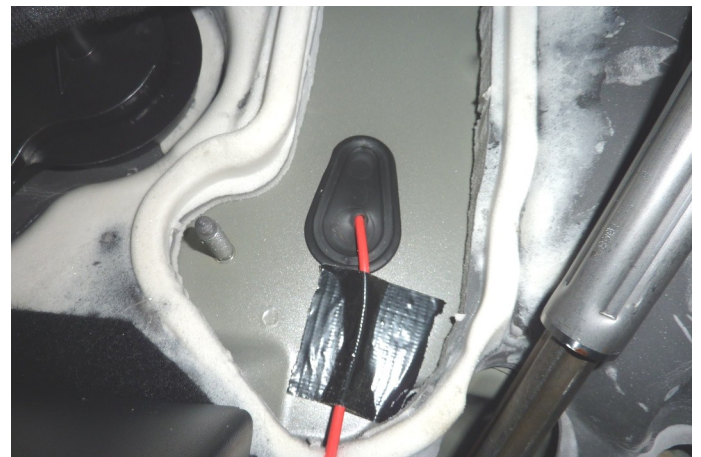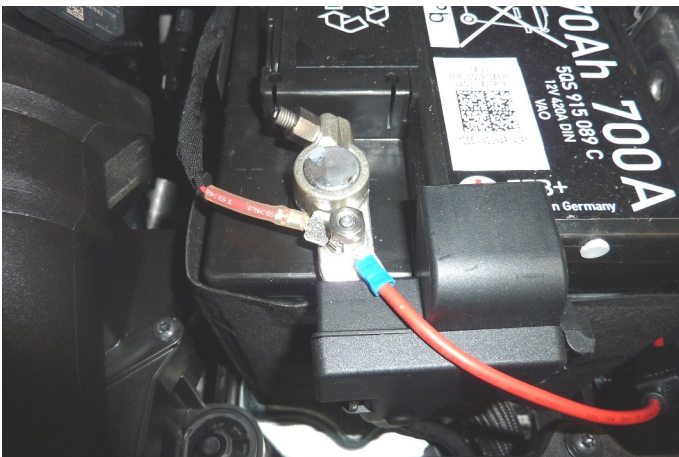


The grommet at the drivers side footwell can be used to route the wire through the bulkhead.



Use an appropriate chassis earth for the black wire and complete the wiring.
Red to Yellow / Black to Green

Fitting Instructions						
Page 8						
Make	Ford	Model	Tourneo Connect	Derivative	MK III 5 Seat	REV B

The fitting instructions show the older style of a-frame. Please use the provided kit to mount the control box, restraint reel and restraint buckle in suitable locations.

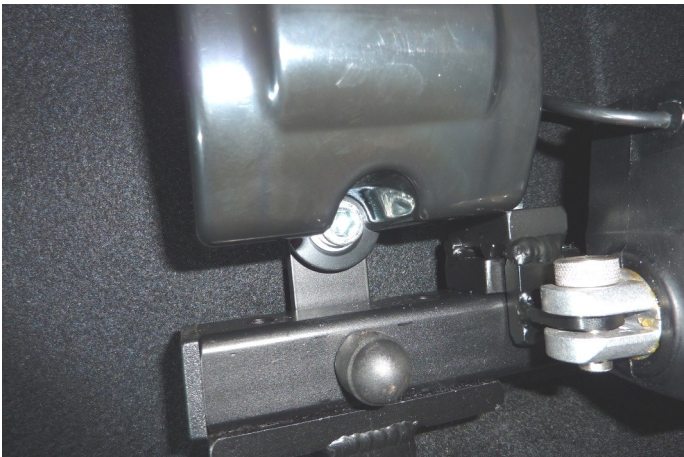

Complete the installation and setup of the hoist.

Fitting Instructions

Page 9

Make	Ford	Model	Tourneo Connect	Derivative	MK III 5 Seat	REV B
------	------	-------	-----------------	------------	---------------	-------



Complete the setup of hoist for the mobility device.

