


Boot Hoist Fitting Instructions

Make	Ford	Model	Transit Custom	Derivative	MK II	REV A
Page 1						




Attach the vehicle specific legs to the a-frame as shown.

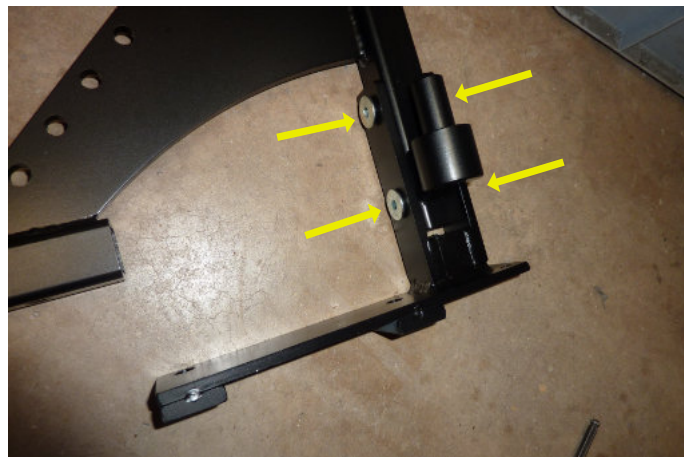

Outer Leg:
4 x M10x20mm Countersunk Set Screws



Upper Leg:
2 x M10x40mm Countersunk Set Screws

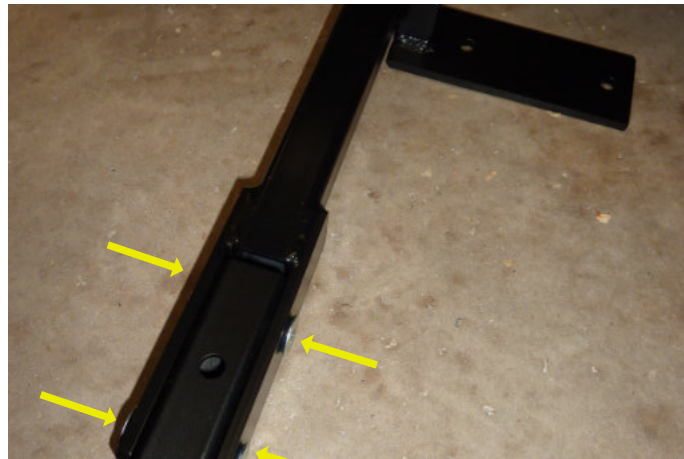

Inner Leg:
4 x M10x20mm Countersunk Set Screws



If ply-lining is fitted, remove the section over the wheel arch.


Boot Hoist Fitting Instructions

Make	Ford	Model	Transit Custom	Derivative	MK II	REV A
------	------	-------	----------------	------------	-------	-------

Page 2



Identify the existing hole through the d-post and open it out to 12mm.

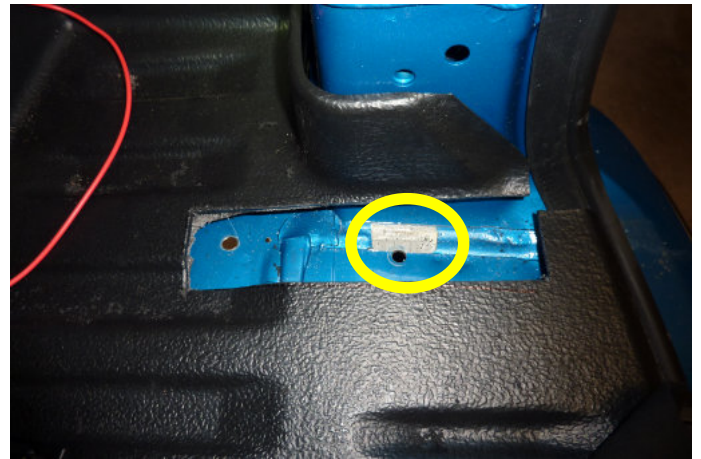

Hold the a-frame so the upper bracket hole is vertically in line with the hole just drilled. The bracket hole will sit slightly higher than the drilled hole due to the flooring.

Make sure the a-frame is level and mark the plastic flooring at the outer leg and then cut it out so the outer leg can sit flush to the vehicle floor.

The small section of floor seal circled will need to be removed to allow the leg to sit flat.



On the design vehicle the 12V battery was in the engine bay.




Run the power cable and attach it to the battery via a fuse holder.

Boot Hoist Fitting Instructions

Make	Ford	Model	Transit Custom	Derivative	MK II	REV A
------	------	-------	----------------	------------	-------	-------


Page 3

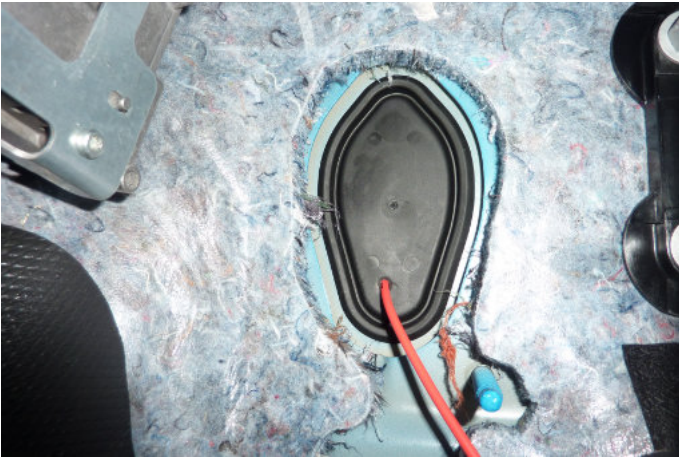

There was a suitable grommet to go through in the bulkhead on the driver side.



A suitable chassis earth was found on the d-post in the design vehicle.

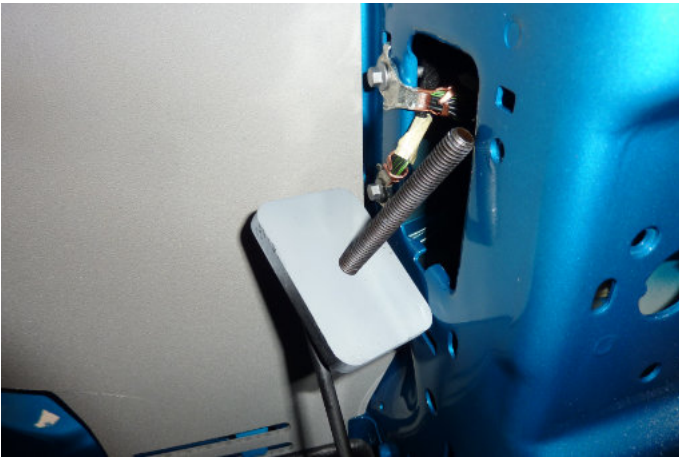

The upper leg is attached through the drilled out hole using the backplate with male thread.



Feed the backplate through the opening as shown to be able to place the thread through the hole from the rear.



Place the spreader plate and pylon over the front of the thread.



Loosely attach the upper leg bracket on the frame to the upper leg bracket on the vehicle. Keep it loose until the other legs are attached to the vehicle.

Boot Hoist Fitting Instructions

Make	Ford	Model	Transit Custom	Derivative	MK II	REV A
------	------	-------	----------------	------------	-------	-------

Page 4



With the a-frame loosely attached to the upper leg, mark and drill the exposed vehicle floor for the 2 fixings shown.



For the rear fixing, place the spreader plate into the cavity shown.



Then feed the male threaded end of the pylon up through the spreader plate and up through the hole previously drilled through the vehicle floor.



Boot Hoist Fitting Instructions

Make	Ford	Model	Transit Custom	Derivative	MK II	REV A
------	------	-------	----------------	------------	-------	-------

Page 5



This is then secured in the vehicle using the M10 nyloc nut and washer.



The inner of the two holes for the outer leg is simply attached using an M10x40mm hex head set screw fed in from above and secured under the vehicle with the spreader plate and nyloc shown.




Ensure an appropriate sealant is applied to any external fixings.

Boot Hoist Fitting Instructions

Make	Ford	Model	Transit Custom	Derivative	MK II	REV A
------	------	-------	----------------	------------	-------	-------

Page 6

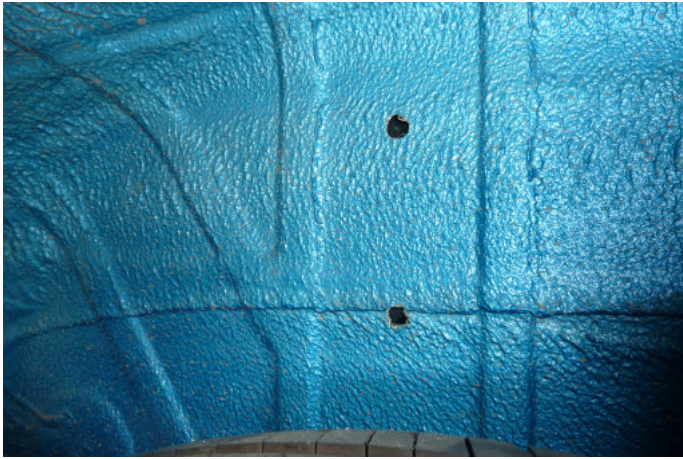

Mark and drill (10mm) the two holes through the wheel arch. These are required for the inner leg.



Apply an appropriate sealant to the inner surface of the plate shown.



Feed the screws through the two holes drilled from the wheel side.



The inner leg then attaches to this plate using two M10 nyloc nuts and washers.

Boot Hoist Fitting Instructions

Make	Ford	Model	Transit Custom	Derivative	MK II	REV A
------	------	-------	----------------	------------	-------	-------

Page 7




This is the view from the wheel arch side externally.



Once in place, the two grub screws of the outer leg can be extended down to contact the vehicle floor for extra stability.



Complete the hoist installation by fitting all parts including the hinge post, arm and actuators. Make sure they are all secure as required.

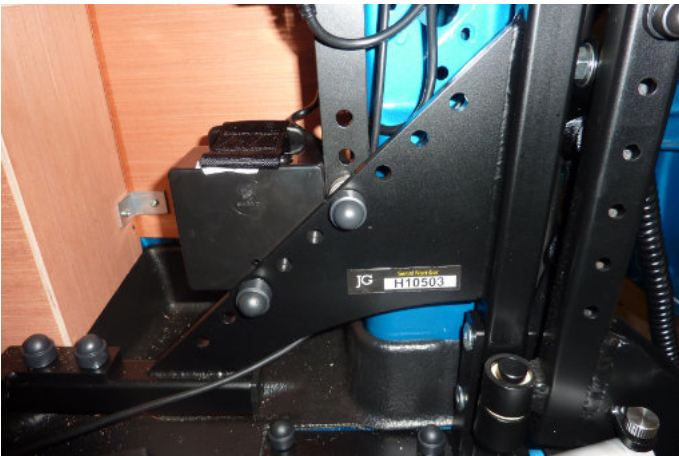

Ensure the control box, seat belt reel and seat belt buckle are mounted in suitable locations. The images show the positions on the design vehicle.

Boot Hoist Fitting Instructions

Make	Ford	Model	Transit Custom	Derivative	MK II	REV A
------	------	-------	----------------	------------	-------	-------



Setup the hoist to suit the mobility device chosen.

