


Boot Hoist Fitting Instructions

Make	Kia	Model	EV3	Derivative	MK I	REV A
Page 1						



Pre-installation pictures.

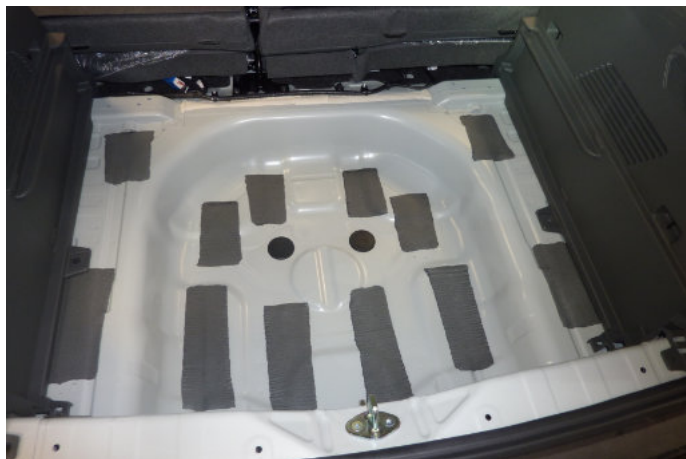

Remove the false floor.



Remove the plastic rivets holding the boot sill in place.



Remove the boot sill.

Boot Hoist Fitting Instructions

Make	Kia	Model	EV3	Derivative	MK I	REV A
------	-----	-------	-----	------------	------	-------

Page 2



Remove these fixings (10mm socket) in order to free the side panel. This can then be flexed out of the way.



Remove the sub-woofer speaker (10mm socket).

Boot Hoist Fitting Instructions

Make	Kia	Model	EV3	Derivative	MK I	REV A
Page 3						



On the design vehicle the 12V battery was in the engine bay.
To access it, unclip and remove the panels shown.



Connect the power to the battery via the fuse holder.



There was a suitable grommet at the drivers side to go through the bulkhead.




A suitable earth was used in the boot as shown.

Boot Hoist Fitting Instructions

Make	Kia	Model	EV3	Derivative	MK I	REV A
------	-----	-------	-----	------------	------	-------

Page 4

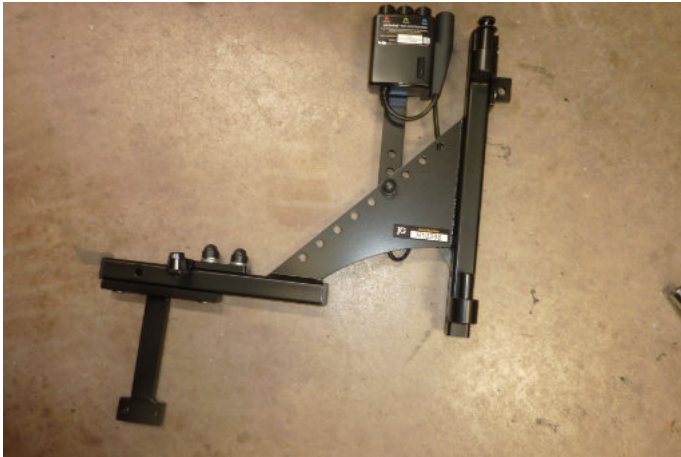

The inner and upper leg portions are attached as shown. It is easier to attach the outer leg to the vehicle first.



Inner Leg:
2 x M10x20mm Countersunk Set Screws



Upper Leg:
2 x M10x 20mm Countersunk Set Screws



Use the blue lines shown as guidance to make a cut out in the side panel for the outer leg.





Boot Hoist Fitting Instructions

Make	Kia	Model	EV3	Derivative	MK I	REV A
------	-----	-------	-----	------------	------	-------

Page 5



Mark the floor with the two lines as shown. Use these lines to determine the position for making the cut as shown for the inner leg.

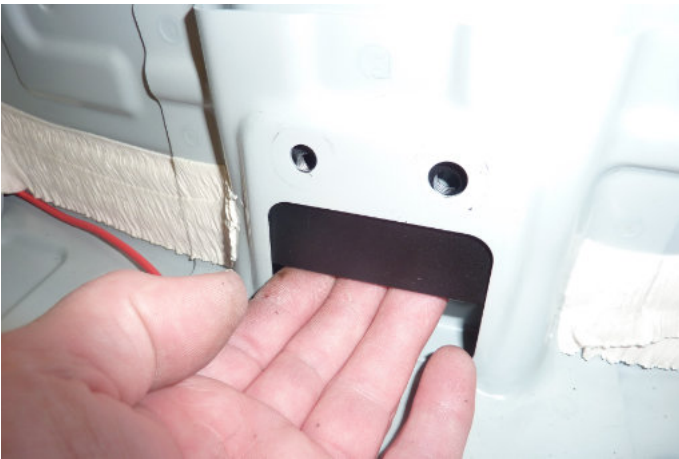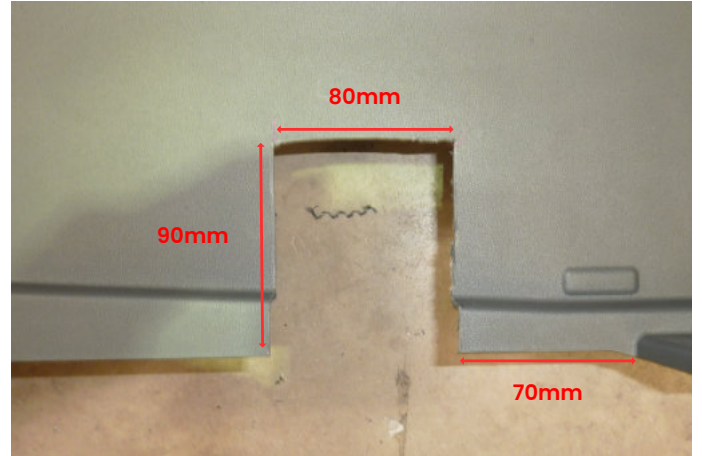

The inner leg back plate can be fed up into the cavity shown and positioned so the holes line up.



Loosely screw on the 50mm pylons ensuring the use of the washers between the pylons and the vehicle.



For the outer leg (not attached to the a-frame yet), slide the backplate into the cavity shown and then attach the leg to the backplate using the M8x25mm Hex Head Set Screws and Washers. Keep these fixings loose for now.



Boot Hoist Fitting Instructions

Make	Kia	Model	EV3	Derivative	MK I	REV A
------	-----	-------	-----	------------	------	-------

Page 6



The upper leg pylon attaches to an existing hole in the vehicle. Drill through the plastic trim at this hole position to allow a hole for the pylon to protrude once attached.



Feed the backplate into the cavity and then screw the upper leg pylon through the hole, into the backplate. Use a spreader plate as shown on the outer side.




Put the plastic trim over the pylon and use a pylon flange to neaten the cut hole.

Boot Hoist Fitting Instructions

Make	Kia	Model	EV3	Derivative	MK I	REV A
------	-----	-------	-----	------------	------	-------

Page 7

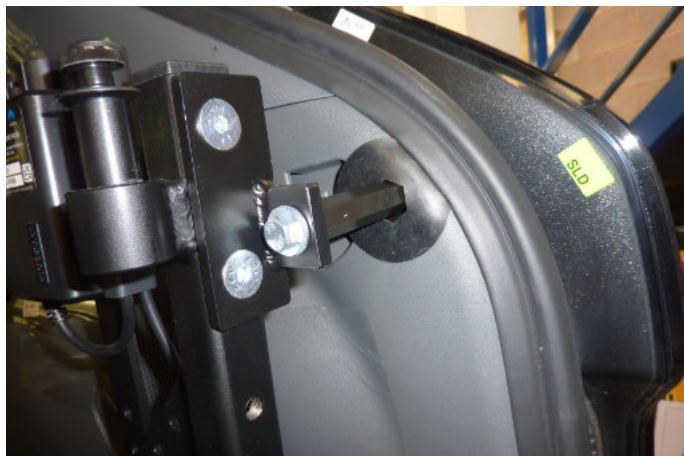

Once the side panel is reattached, the inner and upper legs portions secured to the vehicle can be attached to the a-frame.




Also secure the a-frame in position on the pre-installed outer leg.



All the a-frame fixings can be tightened in situ.

Boot Hoist Fitting Instructions

Make	Kia	Model	EV3	Derivative	MK I	REV A
Page 8						

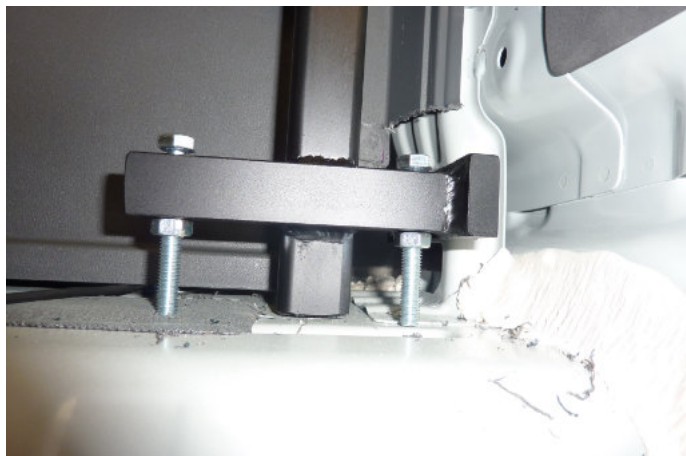


The two M8x60mm hex head set screws can be extended down to the vehicle floor and locked in place to act as stabilisers.



The edge of the boot sill can be trimmed as shown before re-fitting.



Once trimmed, re-fit the boot sill.

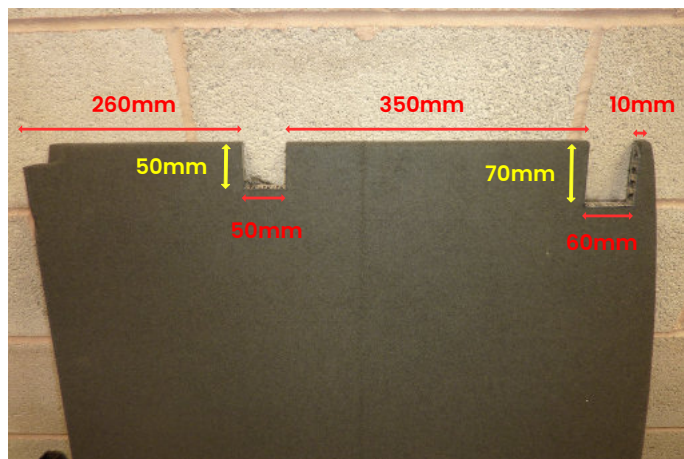

The false floor can be cut using the measurements shown as guidance and then it can be re-fitted.

Boot Hoist Fitting Instructions

Make	Kia	Model	EV3	Derivative	MK I	REV A
------	-----	-------	-----	------------	------	-------


Page 9



Finish the installation by fitting the hinge post, arm and actuators.



Ensure the control box and restraint system is fitting in a suitable location.

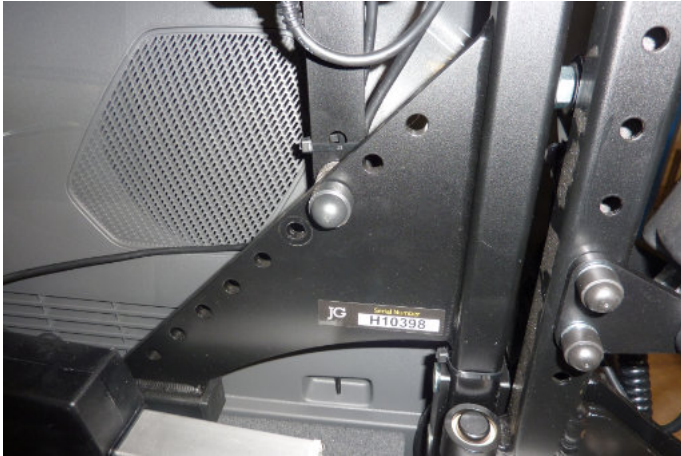

Setup the hoist for the mobility device.



Boot Hoist Fitting Instructions

Make	Kia	Model	EV3	Derivative	MK I	REV A
------	-----	-------	-----	------------	------	-------

