

Fitting Instructions						
Page 1						
Make	Ford	Model	Tourneo Connect	Derivative	MK III Grand 7 Seater	REV A




Fitting Instructions						
Page 2						
Make	Ford	Model	Tourneo Connect	Derivative	MK III Grand 7 Seater	REV A



The 6th and 7th seats must be removed before the hoist is installed. They cannot be used with the hoist fitted.



Remove the two luggage loops using a T40 torx driver.

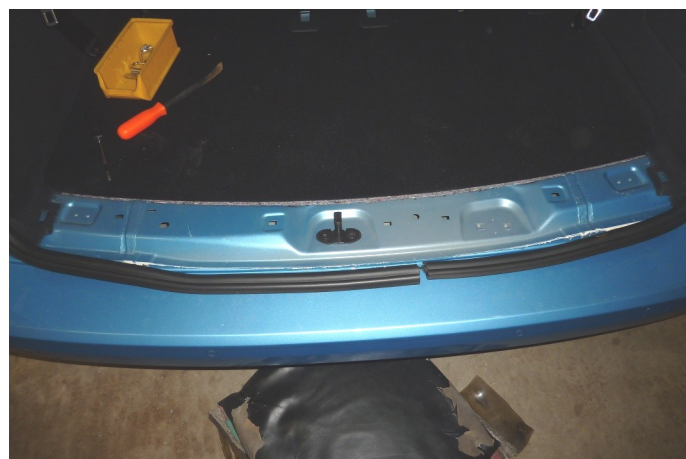

Remove the sill plastic.



Remove the two right hand plastic surrounds at the 7th seat hard points. Here there is a small plastic rivet for which you need a small screwdriver to push down the central pin, then there are two clips holding the other end in.



Partially fold the carpet back.

Fitting Instructions						
Page 3						
Make	Ford	Model	Tourneo Connect	Derivative	MK III Grand 7 Seater	REV A




Remove the screws shown to loosen the side panel.



Run the power wire to the battery which can be found under the bonnet as shown.



A suitable earth point can be found behind the rear side panel.

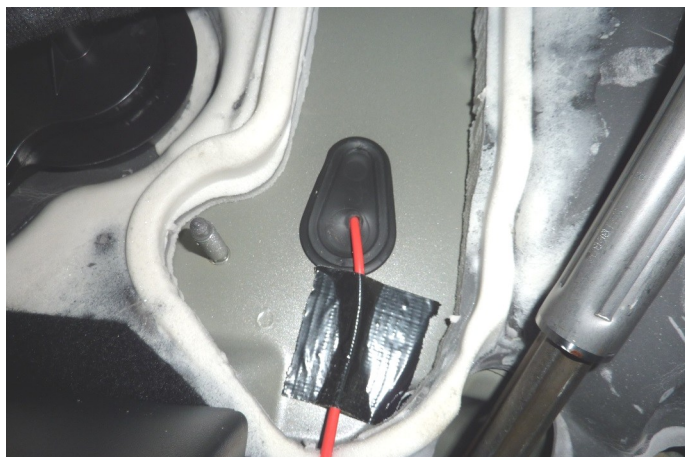

The wire can be routed through the existing grommet in the drivers side footwell.

Fitting Instructions						
Page 4						
Make	Ford	Model	Tourneo Connect	Derivative	MK III Grand 7 Seater	REV A



The outer leg can be used as a template to make the necessary cut out at the right hand side of the sill plastic whilst it is out of the vehicle.




Attach the upper, inner and outer legs to the a-frame as shown.



Inner Leg:
2 x M10x20mm Countersunk Set Screw




Upper Leg:
4 x M10x20mm Countersunk Set Screw

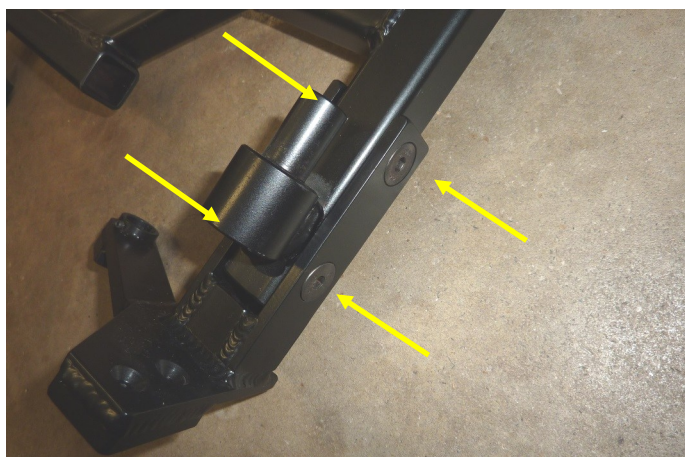

Outer Leg:
4 x M10x20mm Countersunk Set Screw

Fitting Instructions						
Page 5						
Make	Ford	Model	Tourneo Connect	Derivative	MK III Grand 7 Seater	REV A

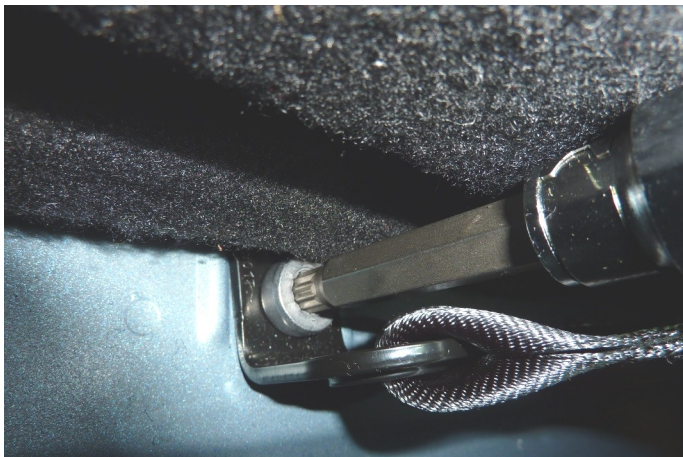

Undo and remove the seatbelt bracket using an M10 star drive bit.



Temporarily attach the a-frame into the vehicle using the 2x M6x30mm countersunk set screws at the outer leg and the 7/16"x40mm countersunk set screw at the inner leg.

Note: the seat belt bracket goes between the floor and the inner leg.



With the a-frame in position, carefully drill a clearance hole through the floor for the outer leg M10 set screw position.

Fitting Instructions						
Page 6						
Make	Ford	Model	Tourneo Connect	Derivative	MK III Grand 7 Seater	REV A



Once the a-frame is positioned, the upper leg location can be established as shown. Drill through the frame to mark the hole position.



Remove the frame and carefully use a hole cutter to open up the trim for the fixing hole.



Note: the side panel only needs to be partially unclipped and flexed sideways slightly in order to gain access to fit the upper bracket.



Attach the upper leg fixing. The threaded backplate can be fed up through the large rectangular opening in the d-pillar.



Use a pylon flange to neatly cover the cut hole.

Fitting Instructions						
Page 7						
Make	Ford	Model	Tourneo Connect	Derivative	MK III Grand 7 Seater	REV A



Re-install the a-frame and attach the inner leg using threadlock.



Attach the outer leg using the two M6x30mm countersunk set screws and M10 set screw. Use the spreader plate under the vehicle with the M10 set screw and apply an appropriate corrosion resistant sealant. The grub screw and ring can be used to assist in holding the hex head set screw in place inside the vehicle whilst tightening underneath.



Fitting Instructions						
Page 8						
Make	Ford	Model	Tourneo Connect	Derivative	MK III Grand 7 Seater	REV A



Attach the upper leg to the protruding pylon.



Make a small cut out in the carpet for the inner leg.



Replace the boot still.



The fitting instructions show the older style of a-frame. Please use the provided kit to mount the control box, restraint reel and restraint buckle in suitable locations.



Fitting Instructions						
Page 9						
Make	Ford	Model	Tourneo Connect	Derivative	MK III Grand 7 Seater	REV A



Complete the setup of the hoist to load the customer's mobility device.

