

Fitting Instructions

Page 1

Make	Nissan	Model	X-Trail	Derivative	MK IV 5 Seater	REV A
-------------	--------	--------------	---------	-------------------	----------------	-------



Fitting Instructions					
Page 2					
Make	Nissan	Model	X-Trail	Derivative	MK IV 5 Seater REV A



Pre-installation



Remove the false floor.



Remove the sub-floor panel.



Remove the boot sill which is held in place with 2 plastic rivets.




Fitting Instructions						
Page 3						
Make	Nissan	Model	X-Trail	Derivative	MK IV 5 Seater	REV A



Remove the 2 M6 screws holding the side panel in place.



Remove the arm rest which is held in place with a single M8 fixing.

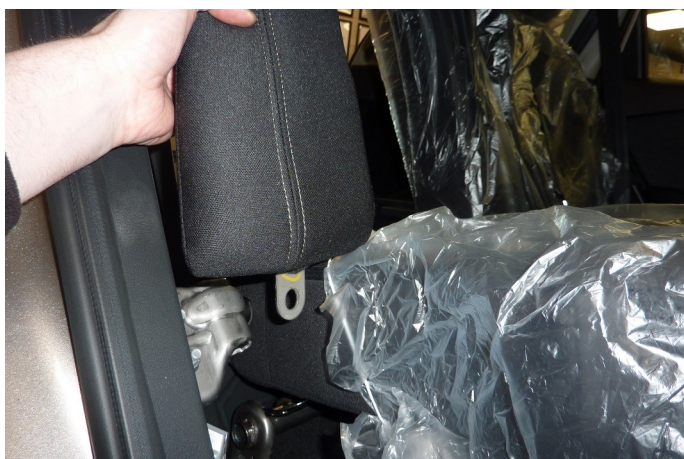

The side panel will then tilt out of the way enough to gain access to the hard points.

Fitting Instructions					
Page 4					
Make	Nissan	Model	X-Trail	Derivative	MK IV 5 Seater REV A



Remove the existing fixings and use the 2x M8x30mm Hex Head Set Screws to attach the inner leg sub-bracket where shown.



This bracket will protrude through the side panel which can be cut to suit. Note: the cut on this image is larger than required.




Cut to the approximate measurements shown.


Fitting Instructions					
Page 5					
Make	Nissan	Model	X-Trail	Derivative	MK IV 5 Seater REV A




Remove the 2 existing screws as shown and attach the outer leg using the 2x M8x45mm Hex Head bolts.

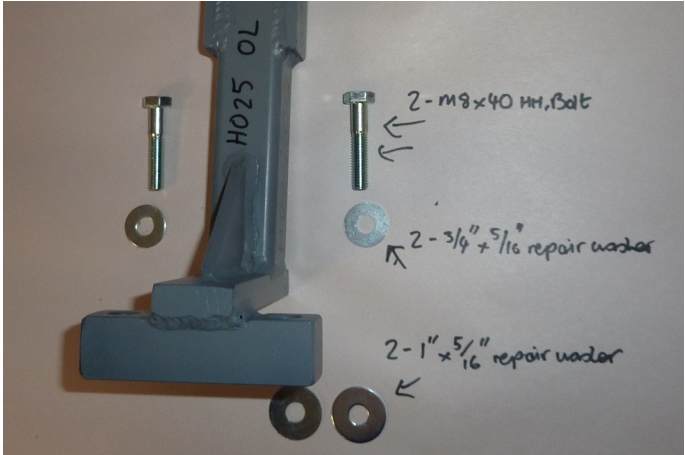


Note the setup of the washers.

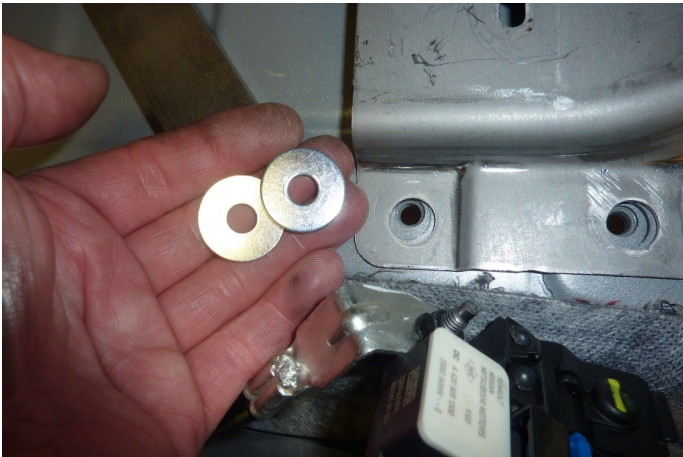

The 2x 1"x5/16" repair washers go under the outer leg where the metal work on the vehicle is stepped down.

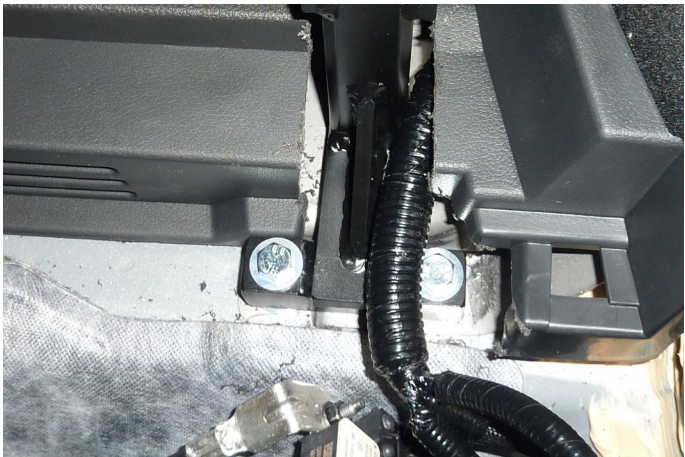

Fitting Instructions					
Page 6					
Make	Nissan	Model	X-Trail	Derivative	MK IV 5 Seater REV A

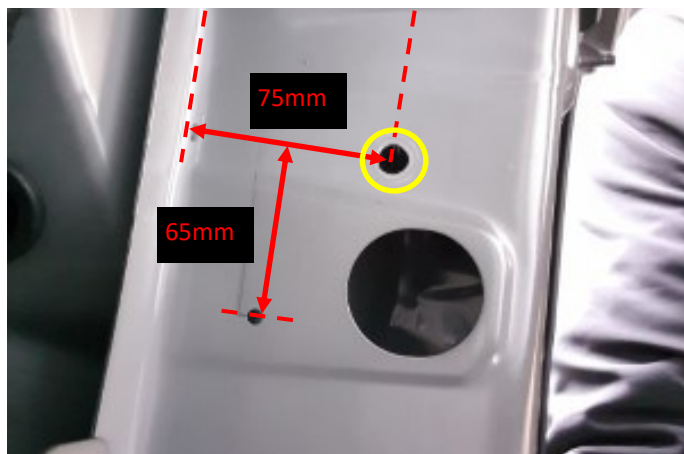

The upper leg attaches as shown.



A 10mm clearance hole is drilled through the metal support and the threaded back plate can be fed through the existing hole shown.

Ensure the use of a drill stop when drilling the side of the vehicle.





The backplate contains an M6x40mm Hex Head Set Screw which can be removed once the plate is tight to the fixings.

Fitting Instructions					
Page 7					
Make	Nissan	Model	X-Trail	Derivative	MK IV 5 Seater REV A




The a-frame can then be attached to the brackets. The upper bracket is shown here.

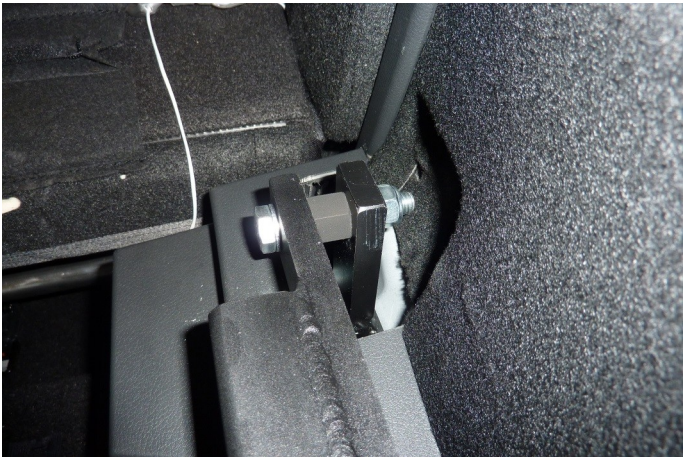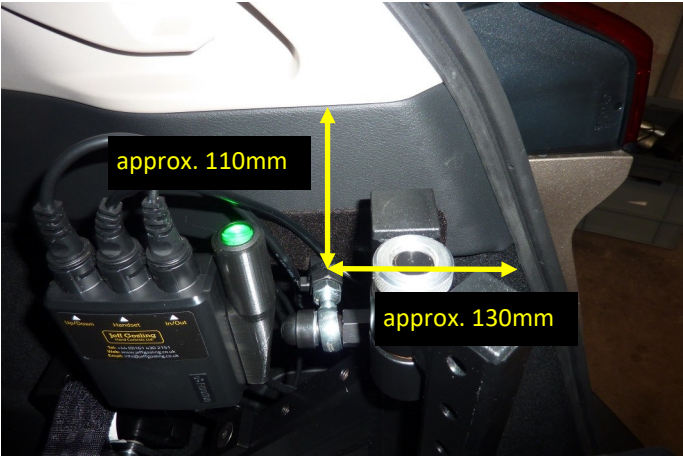

The inner bracket attaches to the sub-bracket via a spacer, washers and M10x60mm Hex Head Set Screw.




The inner bracket attaches to the a-frame using M10x45mm Hex Head Bolts and washers.

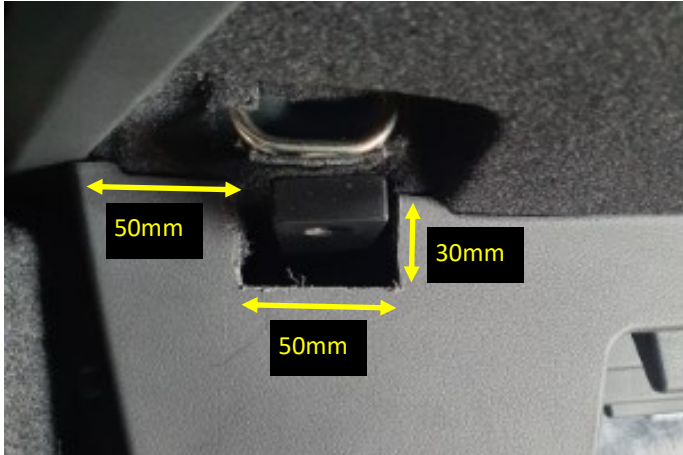


The trim can be cut to accommodate the bracket.

Fitting Instructions						
Page 8						
Make	Nissan	Model	X-Trail	Derivative	MK IV 5 Seater	REV A

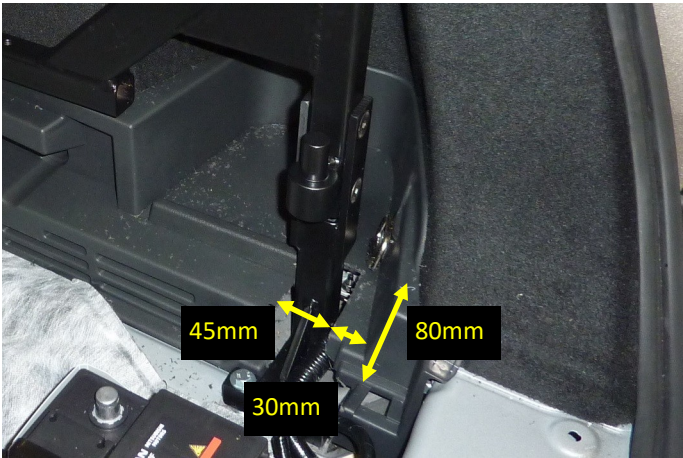

The outer leg attaches via 4x M10x20mm Countersunk Set Screws.
The plastic trim can be cut using the approximate measurements given.

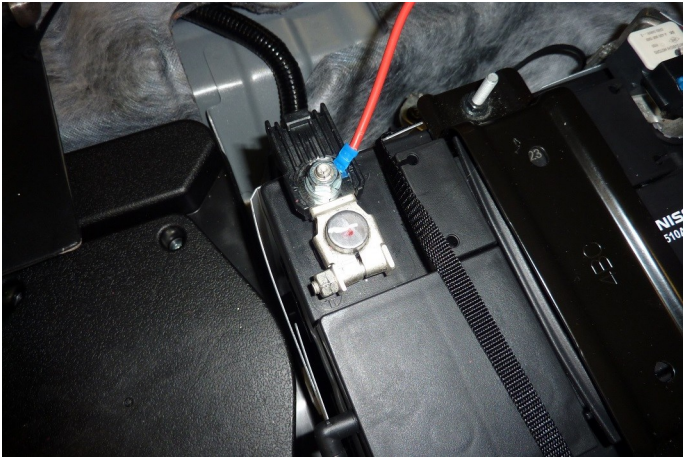

The hoist can be wired to the battery which is located in the boot space.



The power is then attached to the control box as shown.

12V Red to Yellow / Earth Black to Green




A suitable chassis earth should be used for the black wire.



Fitting Instructions					
Page 9					
Make	Nissan	Model	X-Trail	Derivative	MK IV 5 Seater REV A

The fitting instructions show the older style of a-frame. Please use the provided kit to mount the control box, restraint reel and restraint buckle in suitable locations.

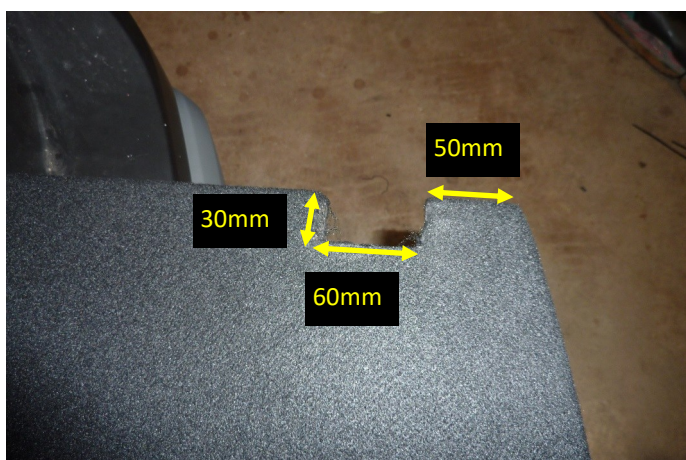

Cut the false floor to fit around the outer leg using the approximate measurements shown.



The false floor can be re-fitted.
Complete the setup of the hoist, ensuring all fixings are secure and it is fully operational.

Fitting Instructions						
Page 10						
Make	Nissan	Model	X-Trail	Derivative	MK IV 5 Seater	REV A



Setup the customers mobility device and the hoist to suit.

